

Tempolat[®]

Temporary composite cement
Dentine-dough

Prescription

Use to cover the drug in the cavity of the tooth.

Composition and characteristics

Tempolat is a white powder temper cement on zinc sulphate(water dentin) base mixed on oil-polymer paste-forming agent. Paste curing in cavity of the tooth upon exposure of moisture during 30-40 min., have good bonding characteristic with ride of cavity and provides hermetic by some kind of broadening in the process of curing. Paste is not toxic, doesn't have the taste, well tolerated by tissue of the tooth and gums. Doesn't contain eugenolum.

Method of application

Paste have to be putted in cavity with the help of burnisher (better slightly heated) over cotton swab 2mm thick coat. Before **Tempolat** will get his hardness recommended to the patient do not eat during one hour. Cured material is easy to take away from the cavity by excavator or probe. The rests of material taking off with passage of bur and

wash with water.

Recommendations


After you take the paste close the cupping glass with lid, checking up availability of consolidate tab.


Conditions of storage and labeling

Condition of storage

Tempolat in packaged form should be stored in closed, cool place, protected from rain and direct sunlight.

Labeling

The  symbol on the package show expire date.

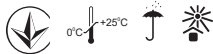

Latus PE, 61030, Textylna str.82, Kharkiv, Ukraine.
Tel +380577332697. www.latus.com.ua